



TruRize™

Clinical Chair



The next step in mobility.

Early patient mobility is a crucial component during the recovery process. Its importance to a patient's health and recovery cannot be emphasized enough.

TruRize helps to promote early patient mobility. It gives caregivers the assurance they need to safely mobilize their patients.

Key Features

- Simple, user-friendly technology helps to equip caregivers with the tools to help promote safe patient handling
- Powered-pivoting design guides the movement of TruRize into positions that help promote patient comfort
- An industry first integrated chair exit alarm system may help to shorten response time to patient falls*



A



B




C



D



E

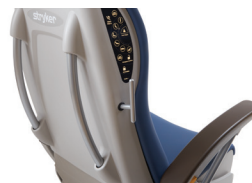

F



- A** Pivoting armrests
- B** Caregiver control panel
- C** Occupant control panel

- D** One-touch brakes
- E** Integrated Foley bag hooks
- F** Vertical transportation handles

Stand Assist Function supports the patient to a standing position, helping to promote early patient mobility.



TruRize™ Clinical Chair



Standard Features

- Powered chair
- Dual-sided Caregiver specific control panels
- Dual-sided lighted occupant control panel with lockout controls
- Six positions (Stand assist, Recline, Upright 2, Upright 1, Transfer, Trend)
- Height Adjustability
- Integrated chair exit alarm system compliant with IEC 60601-1-8 medium priority alarm
- Pivoting armrests
- One-touch brake
- Vertical handles for maneuvering
- Two integrated Foley bag hooks
- Two line management hooks
- 3.5" Casters
- Four fabric color options (Meets CAL 117 flammability)

Optional Features

- Battery Backup
- Nurse Call Cable (1/4" jack)
- Cord reel (domestic only)
- CAL 133 Fire Barrier

Specifications

Model Number	3333-000-000
Overall Width (Upright 1)	34" (86 cm)
Overall Height (Upright 1)	53" to 61" (135 to 155 cm)
Overall Length (Full-Flat)	73" (185 cm)
Seat	
Width	26" (66 cm)
Depth	20.5" (52.1 cm)
Height from floor (Upright 1)	19" to 26" (48 to 66 cm)
Arm Height from Floor	28" to 35" (71 to 89 cm)
Weight Capacity	350 lbs. (158 kg)
Electronics	
Agency Approvals	C-US UL
Volts	100-240 VAC
Ampere Rating	4 A
Frequency	50/60 Hz
Hospital Grade Plug	Yes
Battery (optional)	(2) 12VDC lead acid batteries

Warranty

- One year parts, labor and travel

Stryker reserves the right to change specifications without notice.

Specifications are rounded to the nearest whole number. Conversions are calculated before rounding.

* Powered, clinical chair intended for hospital use only. Compared with the following US or Canada manufactured products: Champion 54 Series, Champion 56 Series, Champion Accent II, Art of Care Powered Bariatric Recliner, Lumex FR588W Series Ger-Chair.

3800 E. Centre Avenue
Portage, MI 49002 USA
t: 269 329 2100
toll free: 800 327 0770
www.patientcare.stryker.com

Mkt Lit-1224 16 NOV 2015 Rev A.1
Copyright © 2015 Stryker
Printed in U.S.A.