


ProForm™


Anatomically placed convoluted foam

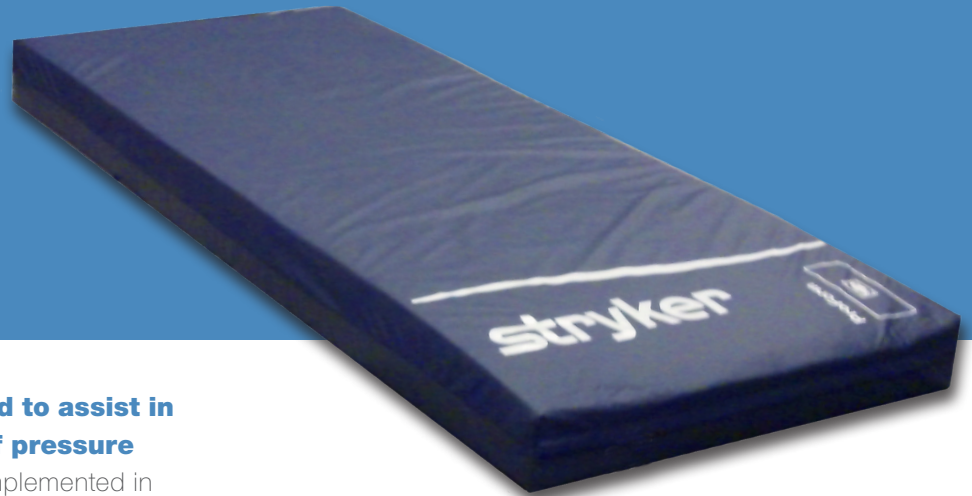


ProForm support surface is used to assist in the prevention and treatment of pressure ulcers and is recommended to be implemented in combination with clinical evaluation of risk factors and skin assessments made by a health care professional¹.

Pressure Redistribution


ProForm's Viscoelastic Foam promotes additional immersion and envelopment enhancing pressure redistribution.

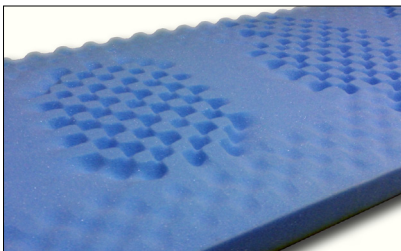

Convoluted foam breaks surface tension at peak pressure points.

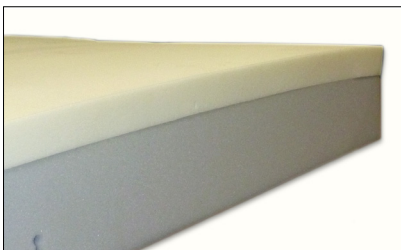

The calcaneus curve supports achilles while keeping calf muscles elongated.

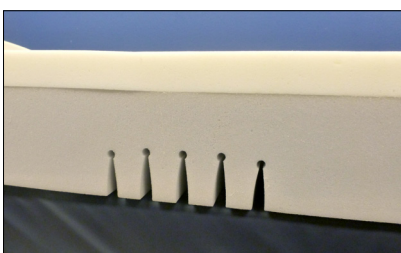

Articulation cuts provide a clean break when HOB is elevated.

ProForm is a non-powered, pressure redistributing support surface through its anatomically placed convoluted foam.

The **anatomically placed convoluted foam** breaks surface tension at peak pressure points.

Supporting areas are non-convoluted to give support to lower peak pressure points.

The **calcaneus curve** utilizes a "Run to the Drop" rather than a distinct angle to support achilles while keeping calf muscles elongated.

Articulation cuts provide a clean break when Head Of Bed is elevated to promote free articulation between support surface and frame.

The **viscoelastic top layer** provides immersion and envelopment to enhance patient comfort and the patient experience.

Covered foam crib enhances patient comfort while maximizing patients surface area.

3800 E. Centre Ave.
 Portage, MI 49002 U.S.A.
 t: 269 329 2100 f: 866-795-2233
 toll free: 800 327 0770

www.stryker.com

ProForm™ Support Surface



Features

- Non-powered design
- Optional fire retardant inner cover meets flammability standard CAL 129, Boston IX-11
- Latex free

Specifications

Model Number	2710
Overall Length	80"/84" (203.2/213.4 cm)
Overall Width	35" (88.9 cm)
Thickness	6" (15.2 cm)
Safe Working Load	500 lb. (226.8 kg)

Options

- Nylon Top Cover
- Dartex Top Cover
- Calcaneus Curve
- Fire Barrier

Warranty

- Three-year, non-prorated warranty on mattress and cover

Stryker reserves the right to change specifications without notice.