

Eye Surgery Stretcher



Enhanced head-end clearance for supreme surgeon access

Ideal for either superior or temporal approach

Positioning controls tuck neatly within the headpiece and allow full articulation

Optional contoured or flat headpiece comfortably supports patient during procedures

Maximum Surgical Access

Stryker's eye surgery stretcher offers the greatest head-end clearance in the industry, giving surgeons ample leg room to maneuver while performing delicate eye procedures. The firmly-locking work surface and mobile design combine to increase patient security and surgical throughput.



Enhanced-clearance, dual-articulating headpiece maximizes surgical access.



Flexible tube doubles as air delivery and drape support.

Accent Color Options



Eye Surgery

Stretcher

stryker[®]



Medical

3800 E. Centre Ave.
Portage, MI 49002
t: 269 329 2100 f: 269 329 2311
toll free: 800 787 9537

www.stryker.com


Standard Features

- 3" (8 cm) Enhanced Comfort Mattress
- 8" (20 cm) Omni Surface™ casters with wheel covers
- Dual cylinder pneumatic backrest
- Dual pedestal hydraulics
- Dual side hydraulic controls
- Fold-down siderails
- Four IV receptacles
- Four-wheel, steel-ring brakes with foot-end brake & steer pedal
- Hands-free Trend.
- Integrated storage tray with oxygen bottle holder
- Retractable fifth-wheel steering
- Rubber bumpers
- Stationary foot-end litter
- Uni-lower pedal
- Dual-articulation, enhanced-clearance headpiece

Specifications

• Model Number	1079	
• Overall Length	90"	(229 cm)
• Overall Width (siderails up)	31.5"	(80 cm)
• Weight Capacity	500 lbs	(228 kg)
• Height Range (to litter top)		
High	36"	(91 cm)
Low	23.25"	(59 cm)
• Litter Positioning		
Backrest	0° - 90°	
Knee Gatch	0° - 30°	
Trend./Reverse Trend.	±18°	
• Headpiece		
Height Adjustment	±4.5"	(11 cm)
Headpiece Tilt	±45°	
• Patient Surface	26" x 87"	(66 x 221 cm)
• Siderails	13" x 55"	(33 x 140 cm)
• Caster Diameter	8"	(20 cm)

Stryker reserves the right to change specifications without notice.

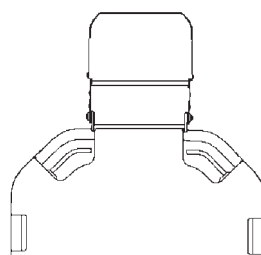

Dual-articulation,
Enhanced-clearance Headpiece

Optional Features

- 3" (8 cm) or 4" (10 cm) Ultra Comfort Mattress
- 4" (10 cm) Enhanced Comfort Mattress
- Pre/post-op shoulder pads assembly (crank backrest only)
- Adjustable arm board and/or clamp (surgery accessory rail required)
- Air delivery system/drape support
- Crank backrest
- Crank knee gatch
- Defibrillator tray
- Defibrillator tray/foot extender
- Temporal or superior wristrest
- Three-sided hydraulic controls
- Footboard/chartholder
- Three-sided brake & steer pedals
- Integrated patient transfer/arm board
- IV caddy
- Restraint straps
- Serving/instrument tray
- Siderail pads
- Surgery accessory rail
- Three-position siderails
- Upright oxygen bottle holder
- Oxygen bottle retaining ring
- Serving/instrument tray holder
- Head pads—Ultra Comfort
 - 3" (8 cm), 2" (5 cm) or 1" (2.5 cm) flat
 - 3" (8 cm) or 2" (5 cm) concave
- IV poles
 - Permanent folding two- or three-stage at foot end
 - Removable
- Extended warranty available

Copyright © 2005 Stryker
LSS-455-03D
PC/SCG-7.5M-1206
Mkt Lit-63-01262005 Rev D